

## Z'-LYTE® Assay Setup Guide on the BMG LABTECH OPTIMA Microplate Reader

NOTE: The BMG LABTECH OPTIMA Microplate Reader was tested for compatibility with Invitrogen's Z'-LYTE® Assay using the Z'-LYTE® Tyr6 kit (PV4122) against JAK2 JH1/JH2 and JAK2 JH1/JH2 V617F kinases. The following document is intended to demonstrate setup of this instrument and provide representative data. For more detailed information and technical support of Invitrogen assays please call 1-800-955-6288, select option "3", then extension 40266. For more detailed information and technical support of BMG LABTECH instruments or software, please contact BMG LABTECH at 1-877-264-5227 or [www.bmglabtech.com](http://www.bmglabtech.com).


### A. Recommended Optics

	wavelength (nm)	BMG LABTECH Filters
Excitation	400 (or similar)	*contact BMG LABTECH
Emission 1	440 (or similar)	*contact BMG LABTECH
Emission 2	520 (or similar)	*contact BMG LABTECH

### B. Instrument Setup


1. Make certain plate reader is turned on, and open up OPTIMA Control software on computer. Insert plate into plate reader.

2. When OPTIMA Control software opens, if you do not have a pre-existing protocol for Z'-LYTE®, select "Test Protocol" from the Test Setup menu bar at the top of the window.

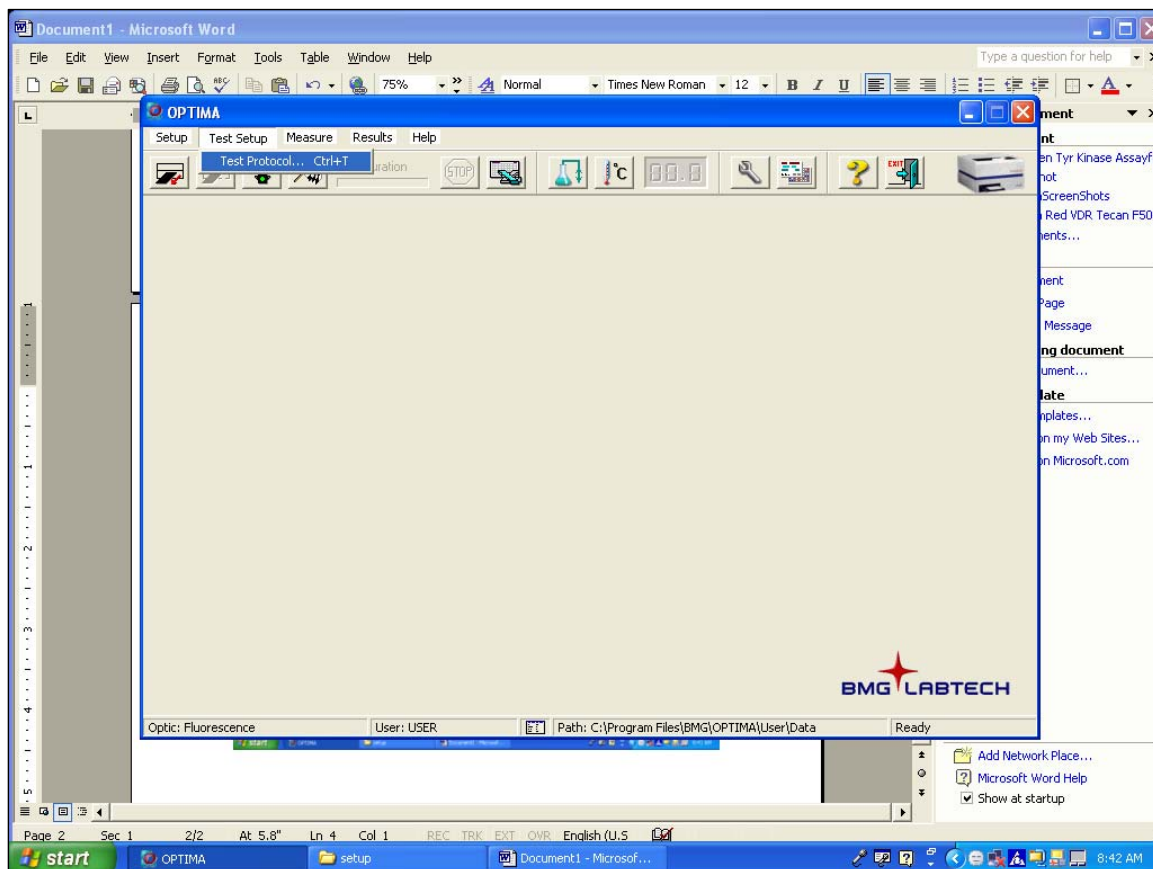


**Setup Guide on the BMGLABTECH OPTIMA Microplate Reader**


3. At this point, a new screen will open (below). Click on the “Show all test protocols” button on the left side of the screen, then select “New” from the tabs at the bottom.

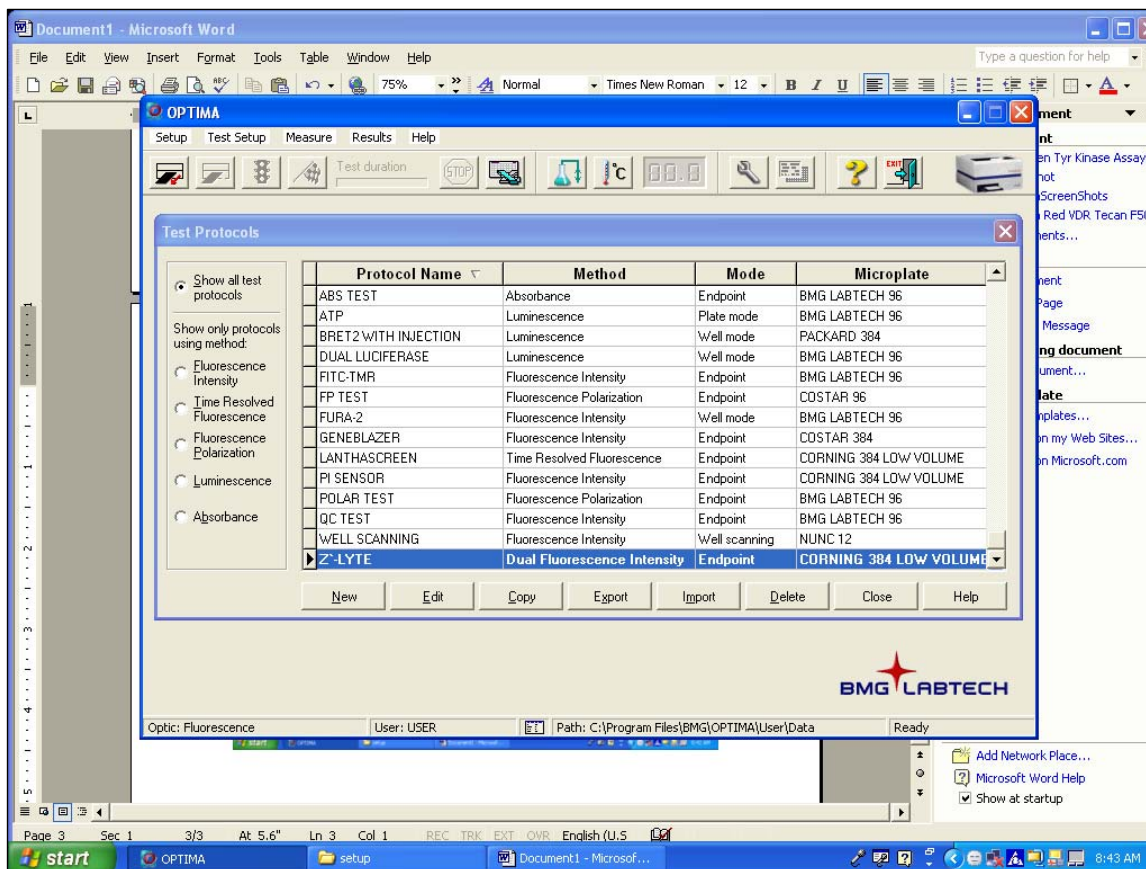


4. A new window will pop up. Select “Fluorescence Intensity” and “Endpoint” and then select “OK”.

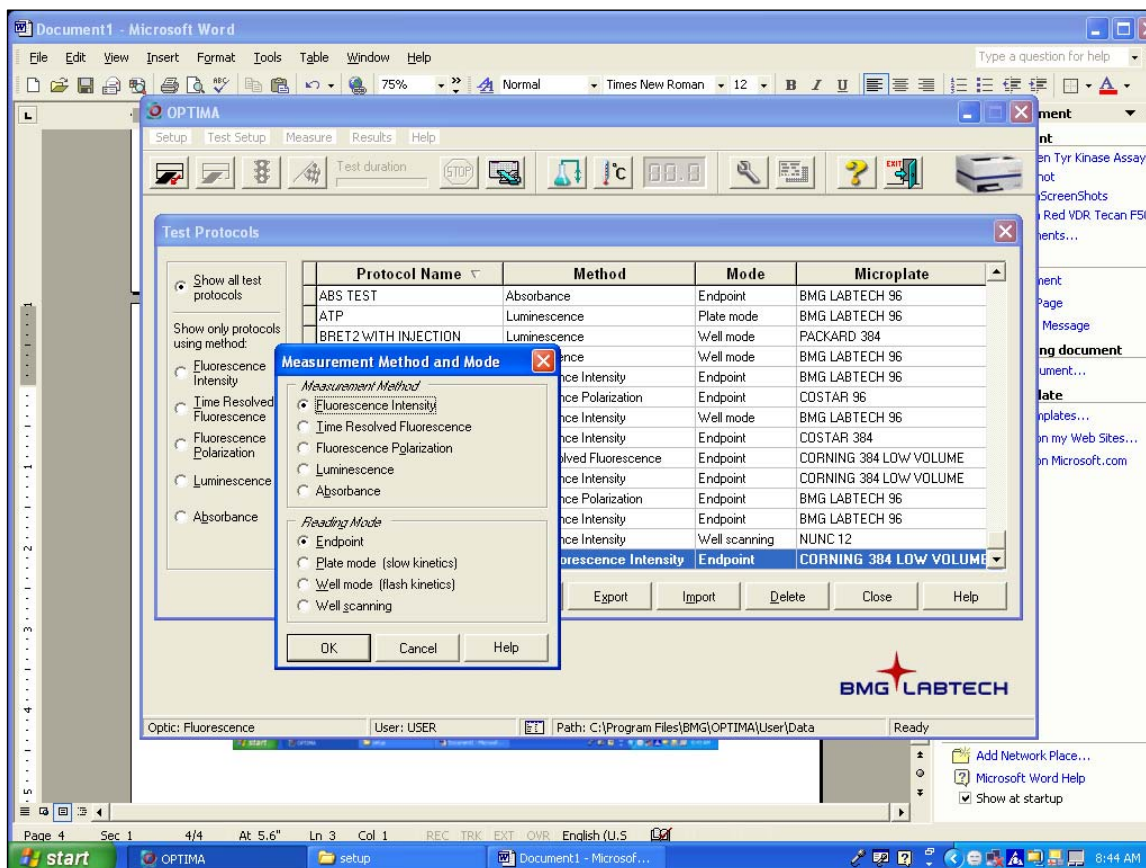

5. A new Protocol window will open automatically. Enter a test name, select plate type, and select the Optic type. Next, set the number of flashes per well and set up the excitation and first emission filter from the drop-down tabs. Finally, under "Filter Settings", enter "2" for the No. of multichromatics. This will allow you to click on the small arrow beside the box.

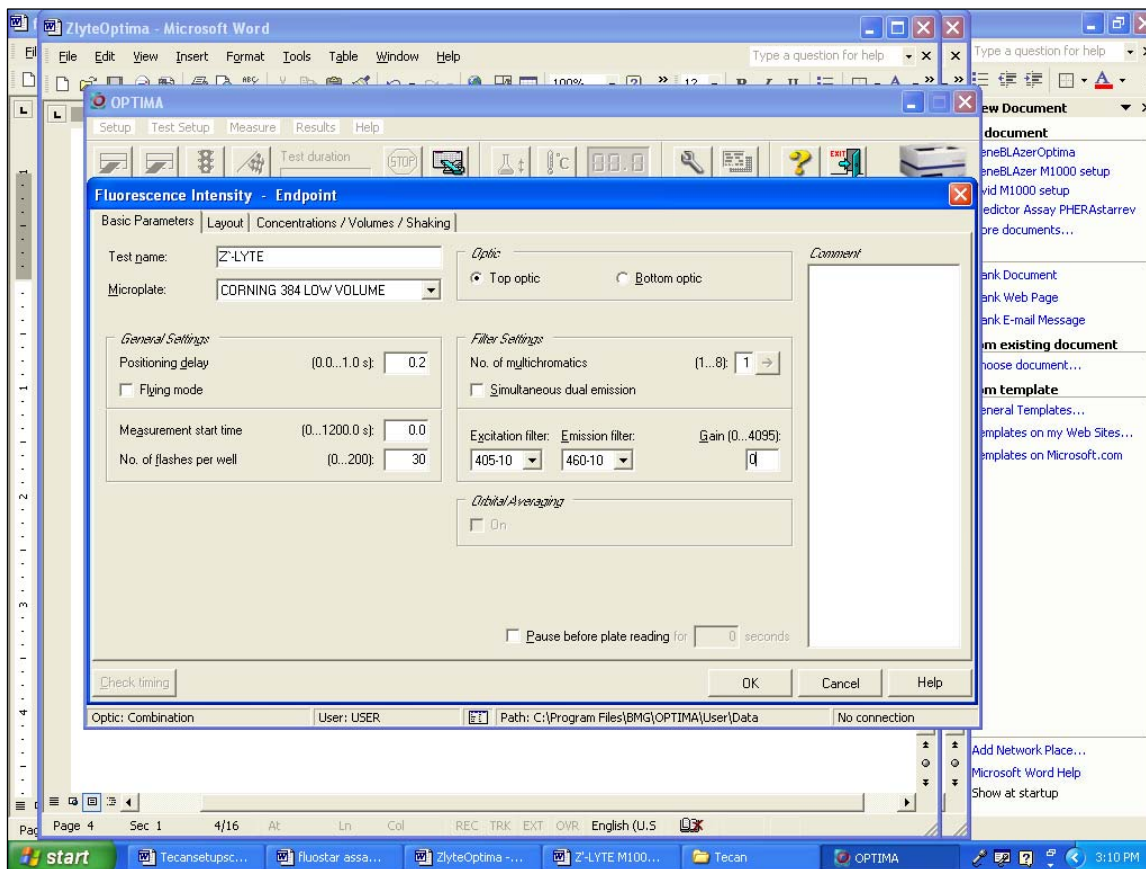


6. Clicking the arrow mentioned in 5 will open another window--select your filters here as shown below, then click "OK" when finished.

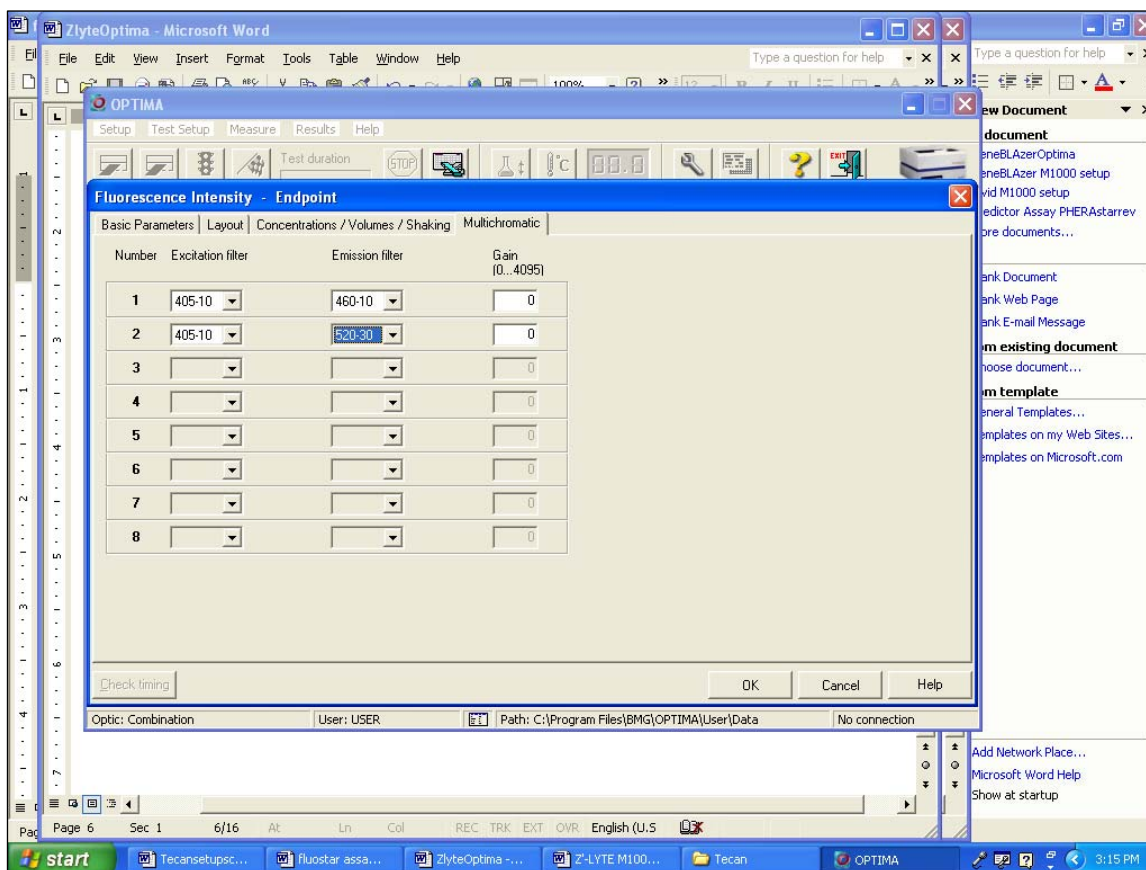




7. You will return to the initial settings window. From the tabs along the top of this window, select "Layout".

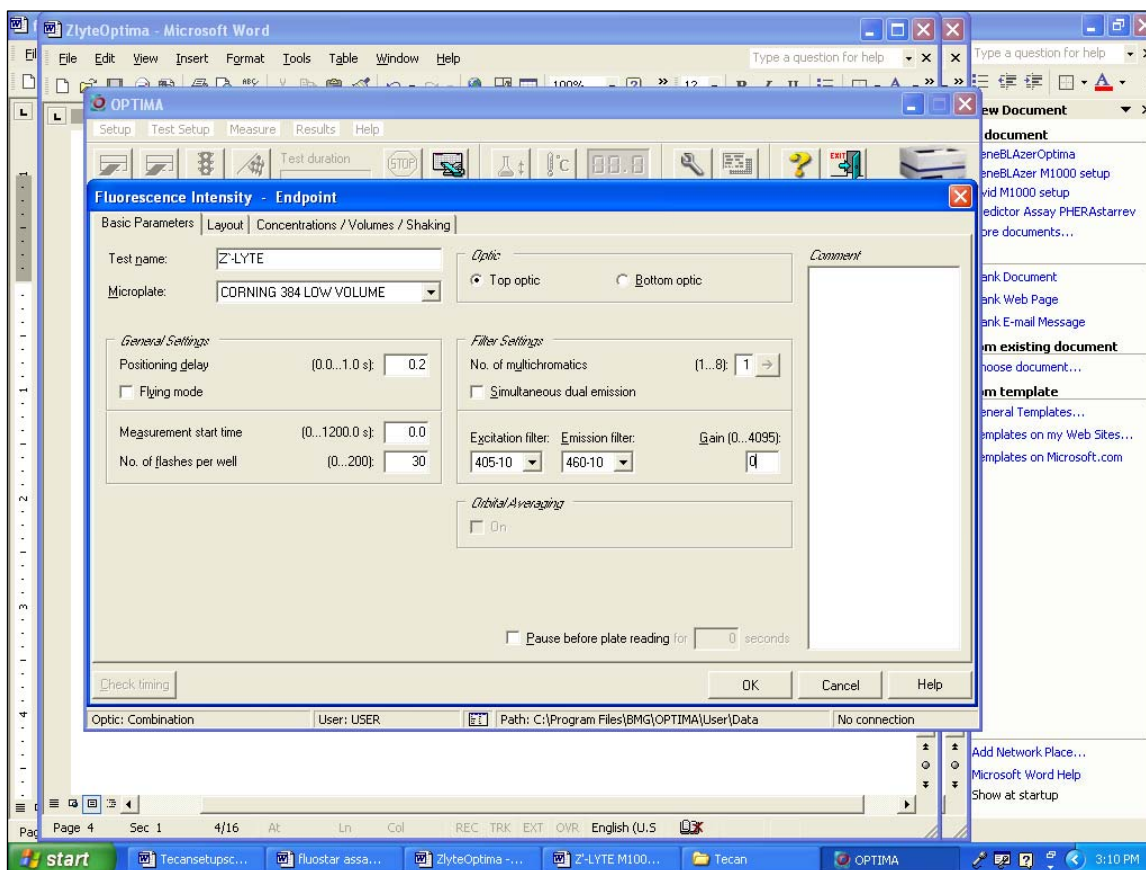

8. A new window will appear. Select the wells you wish to read by highlighting them. When finished, select "OK"

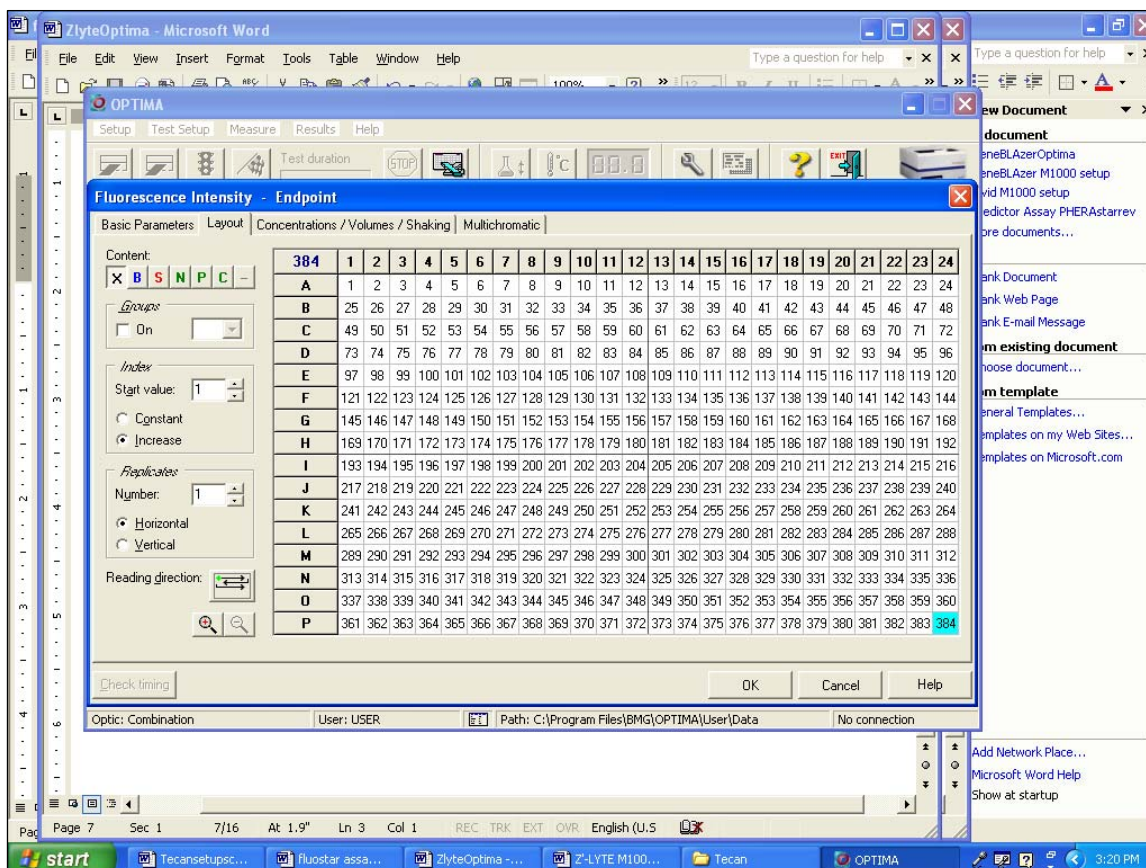



## C. Z'-LYTE<sup>®</sup> Kinase Assay using JAK2 JH1/JH2 and JAK2 JH1/JH2 V617F

NOTE: The following is a sample assay performed for demonstration purposes. The section below describes how the data was obtained, and is not intended for use as an assay protocol. We recommend all first-time users follow the appropriate protocols and/or validation packets provided with their specific assay kits, and include all proper controls. The instrument settings above would be sufficient for any Z'-LYTE<sup>®</sup> assay, the information below is provided as representative data. Assay was run at ATP  $K_{mapparent}$  and a kinase concentration producing approximately 30-40% of maximal phosphorylation, as discussed in Section 9 and 10 of the Z'-LYTE<sup>®</sup> protocols. ATP and kinase concentrations should be optimized for each kinase by the actual user. Specific Z'-LYTE<sup>®</sup> assay protocols and setup information from Invitrogen's own in-house SelectScreen<sup>®</sup> Custom Profiling Z'-LYTE<sup>®</sup>-based kinase assay service can be located at the following link: <http://www.invitrogen.com/content.cfm?pageid=9866>.

1. Prepare initial 100X serial dilution curves in rows A and E of a 384-well plate: Dilute Staurosporine and JAK2 Inhibitor II to a 100X initial concentration in 100% DMSO (100  $\mu$ M). Prepare a set of 1:1 serial dilutions from the initial concentration in a 384-well plate, starting with 80  $\mu$ l in Column 1 and 40  $\mu$ l DMSO in wells 2-20. Add 40  $\mu$ l from well 1 to well 2, and then mix well 2, and take 40  $\mu$ l from well 2 and add to well 3, mix, and so on.

		1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24
Staurosporine	100X	A																							
	4X	B																							
		C																							
		D																							
JAK2 Inh. II		E																							
	4X	F																							
		G																							
		H																							
		I																							
		J																							
		K																							
		L																							
		M																							
		N																							
		O																							
		P																							

**Figure 1: Schematic of initial compound dilution.** Staurosporine and JAK2 Inhibitor II were titrated from a 100  $\mu$ M starting concentration in the initial dilution series by preparing a 1:1 dilution curve in DMSO. A secondary dilution to 4X was then prepared in the rows below the initial dilution curve (lighter gray) using kinase buffer.


2. The 100X serial dilution set is then diluted to a 4X working concentration in Kinase Buffer (50 mM HEPES pH 7.5, 0.01% BRIJ-35, 10 mM MgCl<sub>2</sub>, 1 mM EGTA) in the row below by adding 2  $\mu$ l of diluted inhibitor from the well above to

- 48 µl of kinase buffer. This will produce a final serial dilution starting at 4 µM, which will then produce a final assay concentration starting at 1 µM.
3. Begin to prepare an assay plate: Add 2.5 µL of the compound dilutions per well into a low volume NBS, 384-well plate (Corning Cat. # 3676), in quadruplicate so rows A-D are staurosporine replicates, E-H are JAK2 Inhibitor 2 replicates, etc.
  4. Add 2.5 µl of kinase buffer alone to rows 21 and 22 (0% inhibition no compound control), 23 (0% phosphorylation control, no kinase added) and 24 (Phosphopeptide 100% phosphorylation positive control)
  5. Add 5 µL of the 2X Peptide/Kinase Mixture (2 µM Tyr 06 peptide, 2600 ng/ml JAK2 JH1/JH2 or 1300 ng/ml JAK2 JH1/JH2 V617F, determined experimentally as outlined above) to Columns 1-22. DO NOT ADD TO COLUMN 23 OR 24. Add 5 µL of 2µM substrate alone without kinase to Column 23, rows A-L (0% phosphorylation control) and 5 µl of 2 µM phpsphopeptide control substrate to Column 24, rows A-L (100% phosphorylation control). Add 5 µl kinase buffer alone to the remaining 8 wells (Columns 23 and 24, rows M-P) as a buffer-only reference.
  6. Add 2.5 µL of 4X ATP Solution (200 µM) per well to all Columns to start reaction.
  7. Shake assay plate on a plate shaker for 30 seconds.
  8. Incubate assay plate for 60 minutes at room temperature.
  9. Add 5 µL of the Development Reagent Solution to each well. Use the lot-specific dilutions indicated on your CoA as dilution may vary based upon Z'-LYTE® peptide and Development Reagent A lot.
  10. Shake plate again on a plate shaker for 30 seconds.
  11. Incubate for 60 minutes at room temperature.
  12. Read and analyze as directed in the protocol.

		1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	No Inh.		0	100
Staurosporine	A																								
	B																								
	C																								
	D																								
JAK2 Inh. II	E																								
	F																								
	G																								
	H																								
	I																								
	J																								
	K																								
	L																								
	M																								
	N																								
	O																								
	P																								

**Figure 2: Assay Plate Schematic.** Compound titrations shown in Columns 1-20, Columns 21 and 22 prepared without any inhibitor as kinase activity controls, Column 23 prepared with no kinase (0% phosphorylation) and Column 24 prepared using phosphopeptide control (100% phosphorylation). Note 8 wells in gray in bottom right, which were prepared with out any inhibitor or substrates, as buffer controls.

## D. Results

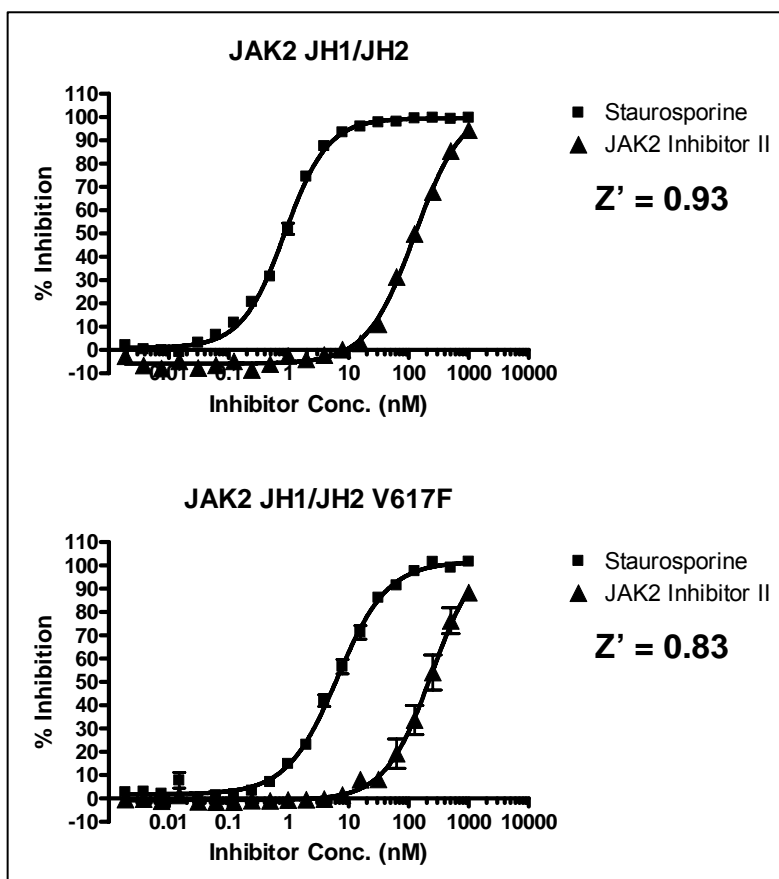

**Figure 1: Z'-LYTE® Kinase Assay.** Z'-LYTE® assay performed using the BMG LABTECH OPTIMA.