

Peroxisomal multifunctional enzyme 2 (MFE2) monoclonal antibody

Cat. no. A21975

Components:	100 µg monoclonal antibody
Lot no.:	See product label
Clone/PAD:	8E5AB11
Isotype:	Mouse IgG3, κ
Gene ID:	3295
Gene Symbol:	HSD17B4
Alternative Names:	Peroxisomal multifunctional enzyme type 2; MFE-2; 17-beta-hydroxysteroid dehydrogenase 4; 17-beta-HSD 4; D-bifunctional protein; DBP; 3-hydroxyacyl-CoA dehydrogenase; 3-alpha,7-alpha, 12-alpha-trihydroxy-5-beta-cholest-24-enoyl-CoA hydratase; DBP; MFE-2; SDR8C1; HSD17B4
Concentration:	1 mg/mL in Hepes-Buffered Saline (HBS) with 0.02% sodium azide as a preservative
mAb PURITY:	Near homogeneity as judged by SDS-PAGE. The antibody was produced <i>in vitro</i> using hybridomas grown in serum-free medium, and then purified by biochemical fractionation.
Reactivity:	Human
Validated Applications:	Immunocytochemistry, Immunoprecipitation, In-Cell ELISA
Suggested Working Concentration:	5 µg/mL for Immunocytochemistry (This is a starting working concentration. The optimal antibody concentration should be determined empirically for each specific application.)
Storage:	Store at 2–8°C. Do not freeze.
Expiration Date:	See product label.

Target Background:

This gene encodes a bifunctional enzyme involved in the peroxisomal beta-oxidation pathway for fatty acids. The gene product also catalyzes the formation of 3-ketoacyl-CoA intermediates from straight-chain and 2-methyl-branched-chain fatty acids. An apparent pseudogene of this gene is present on chromosome 8.



Manufactured exclusively for Invitrogen by MitoSciences, Inc.

For research use only. Not intended for human or animal therapeutic or diagnostic use.


www.invitrogen.com

Life Technologies Corporation • 5791 Van Allen Way • Carlsbad • CA 92008 • Tel: 800.955.6288 • E-mail: techsupport@invitrogen.com

MP21975

MAN0003238

Rev. date: 3 November 2010

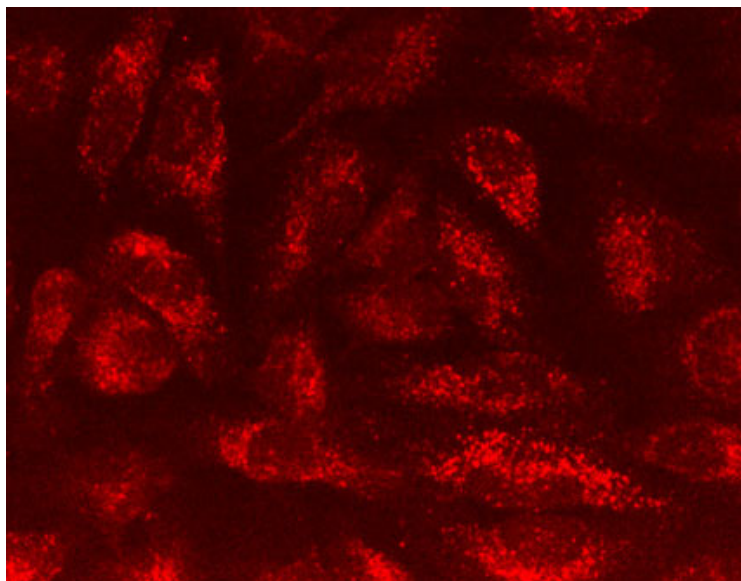

Immunocytochemistry image of Peroxisomal multifunctional enzyme 2 (MFE2) monoclonal antibody.

Human HDFn cells were fixed in 4% paraformaldehyde for 20 minutes and then permeabilized with 0.1% Triton® X-100 for 15 minutes. The cells were incubated with 5 µg/mL of the antibody overnight at 4°C. Alexa Fluor® 594 goat anti-mouse IgG (H+L) was used as a secondary antibody at a 1/1,000 dilution for 1 hour (red). 10% Goat serum was used as the blocking agent for all blocking steps. The target protein localizes mainly in the peroxisome.

The trademarks mentioned herein are the property of Life Technologies Corporation or their respective owners.

Triton® is trademark of Union Carbide.

© Copyright 2010, Life Technologies Corporation. All rights reserved. This information is subject to change without notice.

Important Licensing Information - These products may be covered by one or more Limited Use Label Licenses (see the Invitrogen Catalog or our website, www.invitrogen.com). By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.